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
inFace Sonic Facial Device

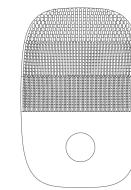Born for love Living for Beauty

inFace Sonic Facial Device is independently developed and manufactured by Shenzhen inFace Technology Co. Ltd., which is affiliated to Shenzhen Migoo Smart Co. Ltd. Since the company was founded, it has been designing the core concepts of technology and intellectual, and striving to explore and innovate

(FEATURES)

Open up a new cleansing frontier of your skin Tleaning — the first sten of skin care


Welcome to use inFace


Elfin Cleanser inFace Sonic Facial Device

Every face wash is a professional enjoyment of deep cleaning gentle care and

1. The centrifugal motor controls the rotating speed (2325pm to 4650pm, and the cushioning design forms the tapping actions at the moment of buffering), and the rotating effects skin directly through sonic vibration technology. Imitating the manual operation helps to absorb drift from deep layer to skin surface, at the same time, avoid thick pores while cleaning.
2. Three-level cleansing massage mode to meet a variety of skin needs.
3. Silicone certified by the FDA (US Food and Drug Administration), and thickened to 4.5mm. The stope design of brush direction, Is degree to 0 degree, fits the factal pore growth direction, which is targeted to clean, no dead angie, well care and do
4. International certified IPXY invalerpoof (not waterproof when charging), it can withstand column watering and could be washed in the water directly.

ACCESSORIES AND

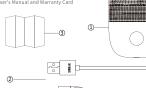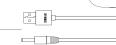
PRODUCT SPECIFICATIONS

Charging Current: 200mA Charging Time: About 2.5h

OPERATION METHOD 1. In the off state, long press the button for 1.5 seconds to turn on the device. The product will automatically shut down if there is no operation within 90 seconds. In the on state, long press the key for 1.5 seconds to turn off the device.

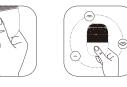3. Touch your face with the brush and move it slowly, stay 2-3 seconds along the Wash off residual foam on your face and continue your daily facial care.

CLEANING METHOD

place in a ventilated dry place.

CHARGING METHOD

2. Shake off the water on the

Notice: Do not use organic solvents, such as alcohol, benzine or thinners, to clean the product.

1. When the battery is low, the white LED above the button will continue to flash.

2. Connect the USB adaptor to the charging interface at the bottom of the device to charge.

3. When charging, the white LED breathing; when full, the white LED is always on and no longer automatically extinguished; unplug the charging cable, the white LED is off. The device cannot be turned on during the whole

PHENOMENON	SOLUTIONS
uring use, the product will utomatically shut down if here is no operation within D seconds.	Software settings, Non-product quality issues, users can turn it on again.
uring use, the irritation caused y sonic vibration may make the kin uncomfortable.	Please adjust the intensity appropriately to use it.
that is the best frequency of sing the product?	Use once a day or every other day(depending on the skin condition to adjust the frequency); if you makeup, it is recommended to use it every day.

(PRECAUTIONS)

(COMMON QUESTIONS)

PHENOMENON	SOLUTIONS
uring use, the product will utomatically shut down if here is no operation within D seconds.	Software settings, Non-product quality issues, users can turn it on again.
uring use, the irritation caused y sonic vibration may make the kin uncomfortable.	Please adjust the intensity appropriately to use it.
hat is the best frequency of sing the product?	Use once a day or every other day(depending on the skin condition to adjust the frequency) if you make up, it is recommended to use it

I Turn of the see year.

I Turn of the power when cleaning the product.

12.In case of danger, do not expose children to the product.

13.Users who had the facial surgery/shaping are recommended to use the product.

13.Users who had the facial surgery/shaping are recommended to use the product.

14. Clean the product in time after use.

15. For the first time to use the product, it is recommended to use up the battery, then fully charge it.

16. The first interpretation right of the instrument belongs to Shenzhen inFace.

If you find any issues during use, please check the abnormality according to the following table. If you cannot solve the problems, contact the after-sales service.

using a non-squeezed foam product, use on-activated device to gently rub eansing cream/powder/glue out of the am, then turn on the device for use.

ventilated place. 2.Keep out of reach of children.

If there is any quality problem after opening the box, please contact the

LIMITATION OF WARRANTY Machines repaired by non-inFace authorized repair center;
 Free gifts and outer packing are not subject to warranty service;

PRESERVATION METHOD

ABOUT WARRANTY

5. Product failure or damage caused by accident or human factors (in improper voltage input, etc.).
If there is any change to the above list, the company's new provisions or th

User Name: Phone Number: Dealer Information Date of Purchase: Order Number: (Fill in the warranty card details when applying for the repairing service. The warranty period is one year for the whole produ

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WARRANTY CARD

AFTER-SALES CONTACT

Service WeChat: migoosmart Working Hours: 9:00-18:00, Monday to Saturday (China time,Day off in Holidays and weekend)

Company: Shenzhen Migoo Smart Co., Ltd. Address: Rm.807, Building U6, Junxiang U8 Intellig Manufacturing Industrial Park, Gushu Community, Xikiang Sub-District, Bao'ar District, Shenzhen, China

X The final interpretation right belongs to Shenzhen Migoo Smart Co., Ltd.